



genetic PCR solutions™

25 years' experience, packed in microtubes

**cfhDNA
MONODOSE
dtec-qPCR Test**

Genetic detection of
Cell-free human DNA

DESCRIPTION

The **cfhDNA MONODOSE dtc-qPCR** tests are individual ready-to-use tubes containing all the components needed to perform a quantitative PCR assay. GPS™ reagents contains BSA and are compatible with all qPCR devices, plate based and glass capillary.

KIT CONTENT

cfhDNA MONODOSE dtc-qPCR tests (INDIVIDUAL TUBES), contains a mixture of specific primers and labelled probe, dNTPs, BSA, polymerase and buffer at optimal concentrations and lyophilized after synthesis. 24, 48 or 96 rxn

Standard Template cfhDNA (RED CAP), $2.4 \cdot 10^7$ target dehydrated copies for positive control.


DNase/RNase free water (GREEN CAP), to dilute standard template. 500 μ l

STORAGE CONDITIONS

All the components of **cfhDNA MONODOSE dtc-qPCR Tests** are stable at room temperature for transport, but should be stored at $-20\text{ }^{\circ}\text{C}$ if not immediately used. Individual tests are stable for one year under this conditions.

PREPARATION OF STANDARD CURVE DILUTION SERIES

Standard Template cfhDNA (RED CAP) contains a high number of template and exist a very significant contamination risk. To minimize the risk of contamination we recommend to open and handle in a separate laboratory environment.

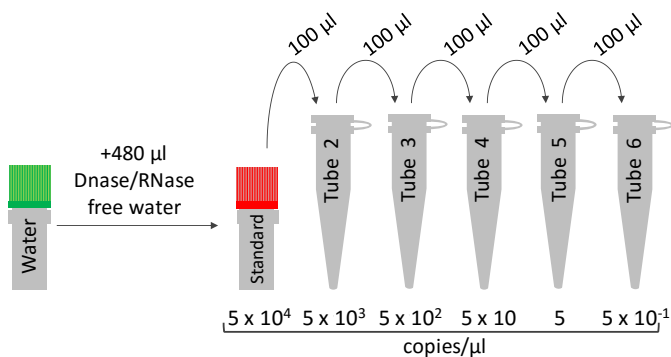

- 1) Pulse-spin the **Standard Template cfhDNA (RED CAP)**, reconstitute with 480 µl of **DNase/RNase free water (GREEN CAP)** and vortex thoroughly, label as number 1
- 2) Pipette 900 µl of **DNase/RNase free water** into five tubes and label as 2 to 6
- 3) Pipette 100 µl of reconstituted **Standard Template**, into tube 2
- 4) Vortex thoroughly and pulse-spin
- 5) Change pipette tip and pipette 100 µl from tube 2 into tube 3
- 6) Vortex thoroughly and pulse-spin
- 7) Repeat steps 5 and 6 with the tubes 4 to 6 to complete the dilution series
- 8) Use 20 µl of each dilution per well to perform the standard curve.

PROTOCOL & AMPLIFICATION REGIME

Add the desired volume of sample ranging from 5 µl up to a maximum qPCR volume of 20 µl and, when needed, complete this final volume by adding DNase/RNase free water (i.e., 7 µl sample + 13 µl water). Vortex thoroughly and pulse-spin. To determine the sample volume, please take into account the possible presence of inhibitors.

IMPORTANT: Protect the mix from prolonged exposure to light.

Cycling parameters are the same for all our MONODOSE dtcc-qPCR assays.

	Step	Time	Temperature
40 Cycles	Activation	60 sec	95 °C
	Denaturation	10 sec	95 °C
	Annealing / Extension and data collection ¹	60 sec	60 °C

¹ Fluorogenic signal should be collected during this step by using the FAM channel

NOTICES, DISCLAIMERS, AND TRADEMARKS

These products included in the present kit are sold exclusively for research use only by the purchaser. They may not be used for human or veterinary in vitro diagnostic (IVD) applications and they may not be re-sold, distributed or re-packaged without express written authorization by GPS™. The developed sequences of primers and probes contained in this mix are a developed industrial property of GPS™. Separate licenses for other than the aforementioned internal R&D applications of the **cfhDNA MONODOSE dtcc-qPCR** may be available. Please inquire via info@genetqpcr.com. PCR is a proprietary technology covered by several US and foreign equivalent patents. These patents are owned by Roche Molecular Systems Inc. and have been sub-licensed by PE Corporation, or licensed by Roche from Life Technologies (formerly Applied Biosystems business group of the Applied Corporation), in certain fields. Depending on your specific application you may need a license from Roche or Life Technologies to practice certain aspects of PCR. In addition, the 5' nuclease assay and certain other homogeneous amplification methods used in connection with the PCR process may be claimed by certain patents of Roche or Life Technologies. Inquiries about obtaining a license under such rights can be made by contacting the Director of Licensing, Life Technologies, 850 Lincoln Centre Drive, Foster City, CA 94404 or the Licensing Department at Roche Molecular Systems Inc., 1145 Atlantic Avenue, Alameda, CA 94501. GPS™ is a trademark of GENETIC ANALYSIS STRATEGIES SL. The purchase of the **cfhDNA MONODOSE dtcc-qPCR** from GPS™ cannot be construed as an authorization or implicit license under any patent rights owned by Roche or Life Technologies.



Genetic PCR Solutions™

by GENETIC ANALYSIS STRATEGIES S.L

CEEI Elche, Pol. Ind. Carrús

Ronda Vall d'Uxó, 125

03206-Elche (Alicante)

☎ Phone: +34-965429901

☎ Fax: +34-966661040

🌐 Web: www.geneticpcr.com

✉ e-mail: info@geneticpcr.com



CENTRO NACIONAL ACREDITADO POR ENAC